



- Extensive Range
- 90+ Years of History
- Extended Aftercare Support
- UKAS Laboratory

Intellicold® Pharmacy Fridges

Free Standing RLDF1524DIGLOCK

A next generation pharmacy refrigerator featuring the latest IntelliCold® 24-Range technology, designed for the safe storage of cold chain medicines in line with Public Health England guidance.

An enhanced IntelliCold® touch screen controller delivers accurate timekeeping, intelligent event logging and Bluetooth connectivity for fast access to temperature data via the Labcold Connect Android App. The flat, easy-clean controller locks automatically to prevent accidental changes, while a factory-fitted digital door lock eliminates the risk of lost or shared keys and puts you in control of authorised access. The user code can be changed at any time, with an optional master key available for added reassurance.

Inside, three independent temperature probes provide advanced protection, including a low-temperature safeguard that prevents freezing. A pre-fitted simulated load has been carefully selected to represent the thermal behaviour of commonly stored pharmacy products. This ensures temperature readings reflect true product conditions, providing meaningful alarms when medicines are genuinely at risk; rather than masking excursions with oversized or unrepresentative loads.

Each unit is supplied with a FREE UK 5-year warranty for complete peace of mind. A UKAS temperature calibration certificate (Lab No. 8898) is available at an extra cost if required.

Specifications

Product Code	RLDF1524DIGLOCK
Gross Capacity	440 Liters
Style	Free Standing
UK Warranty	5 Years
External Dimensions (mm) HxWxD	1875 x 600 x 700
Internal Dimensions (mm) HxWxD	1724 x 480 x 500
Temperature Range	2°C to 8°C
Fridge Set Point	5°C

Specifications (cont.)

Internal Material	Powder Coated Steel
External Material	Powder Coated Steel
Shelves	6
Net Weight	93 Kg
Door Style	Solid
Lock Type	Digital Lock
Digital Temperature Display	Yes
Digital Temperature Control	Yes
High/Low Temperature Alarm	Yes
Min/Max Temperature Recording	Yes
Integral Temp Data Logging	Yes
Automatic Temperature & Fridge Operation Data	Yes
Number of Temperature Probes	3
Internal Light	Yes
Self Closing Door	Yes
Door Open Alarm	Yes
Mains Failure Alarm	Yes
Alarm Battery Backup	Yes
Auto Defrost	Yes
Fan Assisted Cooling	Yes
Remote Alarm Contacts	Yes
Access Porthole	Yes
Access Porthole Size	12mm
Power Consumption Kwh/24hrs	2.55
Ambient Operating Temperature	16°C to 32°C


Features

- FREE UK 5-year warranty
- 24/7 automatic temperature logging which is accessible via Bluetooth (Labcold Connect Android App)
- Fan assisted cooling
- T3 Sensor & Low-temperature cut
- Controller and alarm battery back up
- Power failure alarm
- Minimum/maximum temperature recording and display
- Access port for additional monitoring equipment
- Self-closing door
- Remote alarm contacts
- Fitted with the IntelliCold® touch screen controller
- Event based logging
- Digital door lock with optional master key
- Realistic pre-fitted simulated load
- Door open alarm
- High/Low temperature alarm
- Auto defrost
- Environmentally friendly refrigerant
- Internal light
- Compliant with all leading professional guidelines

Models

	Model	Gross Capacity	Door Style	Lock Type
	Free Standing RLDF1019 Diglock	340 litres	Solid	Digital Lock
	Free Standing RLDG1019 Diglock	340 litres	Glass	Digital Lock
	Free Standing RLDF1519	440 litres	Solid	Key Lock
	Free Standing RLDF1519 Diglock	440 litres	Solid	Digital Lock
	Free Standing RLDF1819	543 litres	Solid	Key Lock
	Free Standing RLDG1819	543 litres	Glass	Key Lock
	Free Standing RLDF1024	340 Liters	Solid	Key Lock
	Free Standing RLDG1024	340 litres	Glass	Key Lock
	Free Standing RLDF1524	440 Liters	Solid	Key Lock
	Free Standing RLDG1524	440 litres	Glass	Key Lock
	Free Standing RLDG1524DIGLOCK	440 litres	Glass	Digital Lock
	Free Standing RLDF1024DIGLOCK	340 Liters	Solid	Key Lock
	Free Standing RLDF1524DIGLOCK	440 Liters	Solid	Digital Lock

Accessories

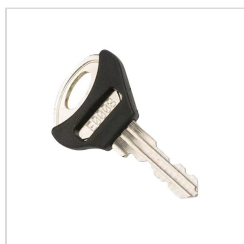

Master Key
LABLOCK3KEYS



Divided Drawer for
RLDF/G1019,
RLDF/G1519,
RLBM1019 and
RLBM1519