



- Extensive Range
- 90+ Years of History
- Extended Aftercare Support
- UKAS Laboratory

Intellicold® Pharmacy Fridges

Undercounter RLDF0524DIGLOCK

A next generation pharmacy refrigerator featuring the latest IntelliCold® 24-Range technology, designed for the safe storage of cold chain medicines in line with Public Health England guidance.

An enhanced IntelliCold® touch screen controller delivers accurate timekeeping, intelligent event logging and Bluetooth connectivity for fast access to temperature data via the Labcold Connect Android App. The flat, easy-clean controller locks automatically to prevent accidental changes, while a factory-fitted digital door lock eliminates the risk of lost or shared keys and puts you in control of authorised access. The user code can be changed at any time, with an optional master key available for added reassurance.

Inside, three independent temperature probes provide advanced protection, including a low-temperature safeguard that prevents freezing. A pre-fitted simulated load has been carefully selected to represent the thermal behaviour of commonly stored pharmacy products. This ensures temperature readings reflect true product conditions, providing meaningful alarms when medicines are genuinely at risk; rather than masking excursions with oversized or unrepresentative loads.

Each unit is supplied with a FREE UK 5-year warranty for complete peace of mind. A UKAS temperature calibration certificate (Lab No. 8898) is available at an extra cost if required.

Specifications

Product Code	RLDF0524DIGLOCK
Gross Capacity	150 Litres
Style	Under Counter
UK Warranty	5 Years
External Dimensions (mm) HxWxD	835 x 600 x 600
Internal Dimensions (mm) HxWxD	690 x 480 x 410
Temperature Range	2°C to 8°C
Fridge Set Point	5°C

Specifications (cont.)

Internal Material	Powder Coated Steel
External Material	Powder Coated Steel
Shelves	3
Net Weight	54 Kg
Door Style	Solid
Lock Type	Digital Lock
Brackets for wall mounting	
Digital Temperature Display	Yes
Digital Temperature Control	Yes
High/Low Temperature Alarm	Yes
Min/Max Temperature Recording	Yes
Integral Temp Data Logging	Yes
Automatic Temperature & Fridge Operation Data	Yes
Number of Temperature Probes	3
Internal Light	Yes
Self Closing Door	Yes
Door Open Alarm	Yes
Mains Failure Alarm	Yes
Alarm Battery Backup	Yes
Auto Defrost	Yes
Fan Assisted Cooling	Yes
Remote Alarm Contacts	Yes
Access Porthole	Yes
Access Porthole Size	12mm
Power Consumption Kwh/24hrs	1.67
Ambient Operating Temperature	16°C to 32°C


Features

- FREE UK 5-year warranty
- 24/7 automatic temperature logging which is accessible via Bluetooth (Labcold Connect Android App)
- Fan assisted cooling
- T3 Sensor & Low-temperature cut
- Controller and alarm battery back up
- Power failure alarm
- Minimum/maximum temperature recording and display
- Access port for additional monitoring equipment
- Self-closing door
- Remote alarm contacts
- Fitted with the IntelliCold® touch screen controller
- Event based logging
- Digital door lock with optional master key
- Realistic pre-fitted simulated load
- Door open alarm
- High/Low temperature alarm
- Auto defrost
- Environmentally friendly refrigerant
- Internal light
- Compliant with all leading professional guidelines

Models

	Model	Gross Capacity	Door Style	Lock Type
	Undercounter RLDG0519 Diglock	150 litres	Glass	Digital Lock
	Undercounter RLDF0524	150 Litres	Solid	Key Lock
	Under Counter RLDG0524	150 Litres	Glass	Key Lock
	Undercounter RLDF0524DIGLOCK	150 Litres	Solid	Digital Lock

Accessories

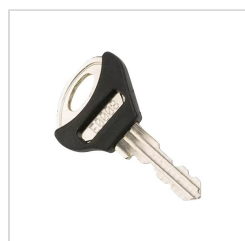

Master Key
LABLOCK3KEYS



Divided Drawer for
RLDF/G0519 and
RLBM0519



Labcold Stacking Kit