


# appleton®




## Low Retention Pipette Tips

# Key features

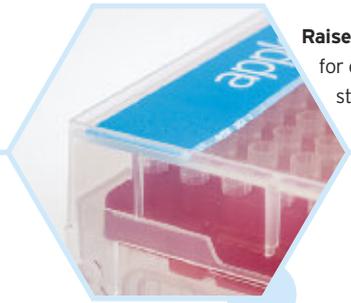


**Raised rim**  
for easy rack  
stacking

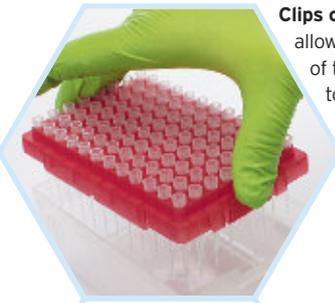


**Clips on side of rack**  
allow safe replacement  
of the rack without  
touching pipette  
tips




**Lot control  
labelling** on  
both racks  
and outer  
boxes

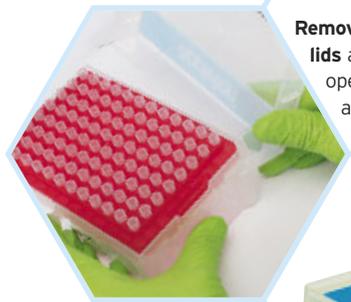

**Removable, hinged  
lids** allow partial  
opening for efficient  
autoclaving

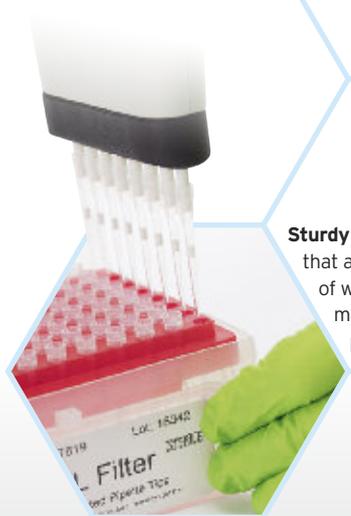

**Sturdy rack and lid**  
that angles out  
of way for easy  
multi-channel  
pipetting

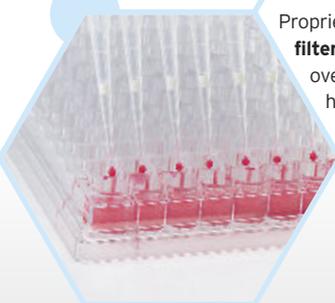


Proprietary **hydrophobic  
filter resin** easily releases  
overdrawn samples and  
has no contaminating  
additives



**Hinged clasp** for  
secure closure  
and easy  
one-hand  
operation

## Low retention tip design

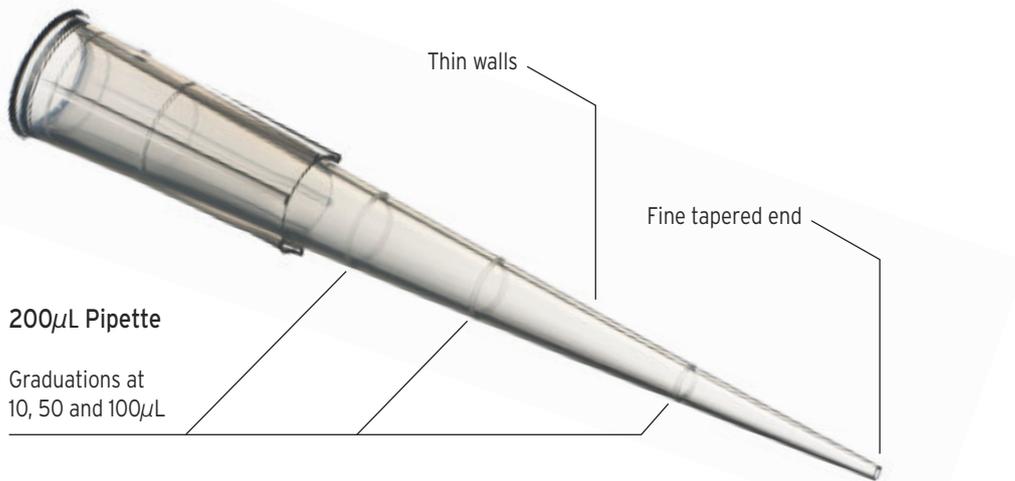

### Low Retention Tips, Appleton

Less fluid retention, less sample loss – Appleton tips prove that less is more when it comes to liquid handling!

Our tips contain a proprietary resin which maximizes low retention. A specially formulated compound increases the hydrophobicity compared to standard tips, enabling viscous fluid in particular to be dispensed completely.

### Clean Liquid Release

To facilitate a cleaner, more efficient liquid release, Appleton tips are designed with thinner walls and a more finely tapered tip construction. This innovative design means less time spent waiting for the final drop to fall off the pipette tip. Users gain complete control working with smaller volumes with no compromise in pipetting accuracy.



## Reloads



1250µL

- Reload decks snap easily into empty tip racks without need for any special equipment
- Tabs on the end of reload inserts are the identical height of the box and do not interfere with multichannel pipettors

- Recessed corners in rack allow steam from autoclaving to easily get into all four corners
- Inserts are imprinted with tip size and lot number data
- Compatible with Appleton standard or filter tip racks



200µL



10µL XL

# Low retention pipette tips

<b>Filter Tips</b>	✓ Low retention ✓ Graduated ✓ Sterile	Description	Code	Pack	Price £
10µL XL		Hinged 10x96 rack	ACT619	960	103.60
20µL		Hinged 10x96 rack	ACT626	960	103.60
100µL		Hinged 10x96 rack	ACT628	960	103.60
200µL		Hinged 10x96 rack	ACT659	960	103.60
1250µL		Hinged 8x96 rack	ACT668	768	81.42

<b>Standard Tips</b>	✓ Low retention ✓ Graduated ✓ Sterile option	Description	Code	Pack	Price £
10µL XL		Bulk bag	ACT610	1000	33.35
		10x 96 tip reload decks	ACT611	960	42.00
		Hinged 10x96 rack	ACT612	960	74.20
		Hinged 10x96 rack, sterile	ACT613	960	58.80
200µL		Bulk bag	ACT621	1000	32.20
		10x 96 tip reload decks	ACT624	960	42.00
		Hinged 10x96 rack	ACT622	960	51.80
		Hinged 10x96 rack, sterile	ACT623	960	57.96
1250µL		Bulk bag	ACT661	1000	30.80
		8x 96 tip reload decks	ACT664	768	42.00
		Hinged 8x96 rack	ACT662	768	51.80
		Hinged 8x96 rack, sterile	ACT663	768	56.00

XL = extended length

For pipettor compatibility please check the website

## Sales & Customer Service

**Tel:** +44 (0)121 458 7740 **Fax:** +44 (0)121 458 5510 **Email:** sales@appletonwoods.co.uk

## Distributed by

**Appleton Woods Limited**

New Lindon House, Catesby Park, Kings Norton, Birmingham B38 8SE, United Kingdom  
www.appletonwoods.co.uk

E&OE



For the Life Scientist