


SORRY

No Image

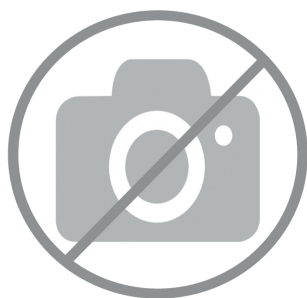

Available

The Labcold Sparkfree range is engineered to store flammable materials safely with no sources of ignition inside the chambers. Anything that could cause a spark is located outside of the chamber.

Warning: Do not store open containers of volatile substances in this freezer.

SORRY

No Image



Available



For extra security all sparkfree models are fitted with a built in lock and digital temperature display. The display features high and low audio and visual temperature alarms as standard in addition to a door open alarm.

RLVF1517 Sparkfree Freezer

Part Number	RLVF1517
Capacity (Litres)	406
Style	Free standing
Exterior Dimensions (HxWxD)	1875 x 600 x 700mm
Interior Dimensions (HxWxD)	1690 x 440 x 500mm
Net Weight	97kg
Energy Consumption	3.16kWh/24hr
Operating Temperature	-18°C to -25°C
Factory Set Point	-20°C
Interior Construction	Painted Aluminium
Exterior Construction	Painted Steel
Door Style	Solid
Door Hinging	Right (Can be reversed at time of manufacture)
Door Lock	✓
Door Open Alarm	✓
Adjustable Shelves	6
Shelf dimensions (WxD)	413mm x 438mm
Fan Circulated	✓
Rollers	4
Temperature Display	✓
High Temperature Alarm	✓
Low Temperature Alarm	✓
Access Port	✓
Automatic Defrost	✗
2 Year Warranty (UK only)	✓